



Dog Articular Chondrocytes (AC)


Catalogue No.: abx700013

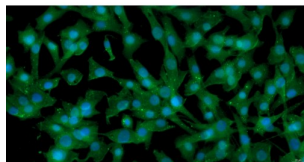

Immunofluorescence identification of collagen II.

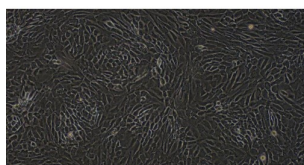

Morphology of P1 Canine Articular Chondrocytes.

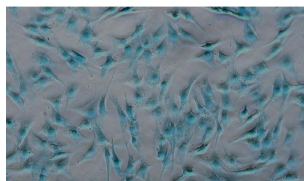

Morphology of Canine Articular Chondrocytes (Alcian blue staining).

Dog Articular Chondrocytes (AC) are Adherent Dog Chondrocytes from Dog Articular cartilage.

Target: Articular Chondrocytes (AC)

Origin: Dog

Host: Dog

Purity: Negative for HIV-1, HBV, HCV, mycoplasma, bacteria, yeast and fungi.

Storage: Shipped at -70 °C. Upon receipt, store in liquid nitrogen (-196 °C). Avoid repeated freeze/thaw cycles.

Validity: 12 months.

Buffer: Contains 90% FBS and 10% DMSO.

Biological Activity: Cell activity: > 85% (viability by Trypan Blue exclusion)

Note: THIS PRODUCT IS FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC, THERAPEUTIC OR COSMETIC PROCEDURES. NOT FOR HUMAN OR ANIMAL CONSUMPTION.

This product is shipped with dry ice.

Directions for use: Recommended Cell Culture Conditions:

DMEM + 10% FBS + Dog Chondrocyte growth supplement + 1% Penicillin-Streptomycin Solution @ 37 °C, 95% air, 5% CO₂.

Cell Recovery: Thaw cells in a 37 °C water bath with shaking until the mixture has dissolved. Transfer to a centrifuge tube and add culture medium (see Recommended Cell Culture Conditions above) at a volume 3-5 times the volume of the cells. Centrifuge at 1000 RPM for 5 minutes and discard the supernatant. Transfer to a T25 flask for culture.

Suggested Cell Passage Procedure: Cells should be 85-95% confluent before cell passage is carried out.

1. Discard the medium and wash with PBS 1-2 times.
2. Add 1 ml of Trypsin at 37 °C, then observe the cells under a microscope.
3. When the cells appear retracted and rounded, gently tap the culture flask to detach the cells. Stop the trypsinization by adding 2 ml of culture medium containing 10% serum.
4. Add fresh medium to resuspend the cells. The recommended ratio of primary cells is 1/2. Pipette to obtain a single cell suspension.